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TABLE OF CONTENTS

▪ INTRODUCTION	_____	p. 23
▪ SAFETY PRECAUTIONS	_____	p. 23~24
▪ PRECAUTIONS BEFORE USE	_____	p. 25
▪ PRECAUTIONS IN USE	_____	p. 26
▪ PRECAUTIONS OF MAINTENANCE AND INSPECTION	_____	p. 27
▪ NAME OF PARTS	_____	p. 28
▪ EXPLANATION OF FUNCTIONS	_____	p. 29
▪ HOW TO USE		
1. Basic Operations	_____	p. 30
2. Connect to dedicated app	—	p. 31
3. Setting	_____	p. 32~33
4. LED	_____	p. 34
5. Treatment	_____	p. 35
▪ TREATMENT PROCEDURE		
1. Basic Procedure	_____	p. 36
2. Face	_____	p. 37~38
3. Body	_____	p. 39~40
▪ PRODUCT SPECIFICATIONS	_____	p. 41
▪ ASSURANCE PROVISION	_____	p. 42

INTRODUCTION

- Thank you for purchasing Zeus III.

In order to use the device properly and safely, it's necessary to read this instruction manual carefully and follow all safety precautions. Please include this instruction manual if you borrow or give to another person.

- This instruction manual contains the warranty card.

Confirm the information such as "The date of purchase" and "The store name or dealer name", and store it carefully.

SAFETY PRECAUTIONS

24

Be sure to read this instruction manual carefully before using, and always use properly.

To avoid the risk of injury, loss of life, and damage to any property, the following symbols are used. The marks and meanings of symbols are shown as below.

 DANGER	DANGER indicates an imminently hazardous situation which, if not avoided, could result in death or serious injury.
 WARNING	WARNING indicates a potentially hazardous situation which, if not avoided, may result in death or serious injury.
 CAUTION	CAUTION indicates a potentially hazardous situation which, if not avoided, may result in minor injury or property damage. ※Property damage means damage to household belongings, livestock or pets.

■ MEANINGS OF SYMBOLS

	This symbol indicates prohibition. ※It is indicated a special instruction inside or near a picture.
	This symbol indicates compulsory actions.
	This symbol indicates caution.

SAFETY PRECAUTIONS

WARNING

To avoid accidents and discomfort.

Prohibition

■ Do not use this product with any of the following medical electrical equipment.

- Malfunction of medical electronic equipment causes accidents or discomfort.
- Implantable medical electrical equipment such as a pacemaker.
- Life support medical electrical equipment such as an artificial heart-lung device.
- Wearable medical electrical equipment such as an electrocardiograph or a dropper.

■ The following persons must not use the product.

- Persons with heart disease, abnormal blood pressure or malignant tumors.
- Persons with diseases with fever, infectious diseases or acute diseases.
- Persons who are receiving dental treatment.
- Persons with weakened immune systems or who use immune suppressants due to immunity disorders like AIDS etc.
- Persons with hemophilia, bleeding disorders, diabetes or facial neuralgia.
- Persons with paresthesia or dysalgnesia.
- Persons with Keloid diathesis, facial melanosis, photosensitivity or pigmentation disorders.
- Persons with capillary dilatation due to long-term steroidal hormone treatments or liver dysfunction.
- Persons who have allergies, atopic dermatitis, skin diseases, sensitive skin or weak skin.
- Persons who have epileptic stroke or discomfort due to flashing lights.
- Persons who have aftereffects of head or ophthalmologic aftereffects due to traffic accidents.
- Persons who have had eye operations (in the case of LASIK, consult with the doctor before use).
- Persons who are pregnant, on period or breast-feeding.
- Children or the person who cannot express their feelings in words.
- Persons under the influence of alcohol or with a hangover.
- Persons who take drugs that cause drowsiness.
- Persons who are extremely tired or do not feel well.
- Persons who have had skin problems such as spots or rashes due to artificial fibers.
- Persons who have or have had skin-inflammation due to cosmetics.
- Persons immediately after having chemical peeling at beauty salons etc.
- Persons who have had skin problems after using beauty equipment.

25

Prohibition

■ Do not use on the following parts.

- Except specific parts.
- Around heart, eyeballs, upper eyelids, lips, and mucous membrane such as inside the mouth.
- Private parts, anus, around the Adam's apple, and ears.
- Parts with morbid spots, moles, and viral warts.
- Tingled and inflamed parts due to excessive sunburn.
- Parts with wounds, eczema, psoriasis, herpes, purulent acnes and/or pimples.
- Parts with pathological spots due to itching, flushing, physical stimulus.
- Parts where surgery was performed.
- Any parts immediately after hair removal or shaving has been performed.

Prohibition

※If you are seeing a doctor, consult with the doctor before use.

PRECAUTIONS BEFORE USE

WARNING

26

Prohibition

- Do not use this product with other beauty equipment at the same time.
- Do not use this product for any purpose other than beauty and under conditions of malfunction or damage. This product is a treatment equipment specialized for use in beauty. Do not use for other purposes such as therapeutic purposes.
- Do not look at or directly view the LED light.
- Do not use continuously on the same parts. Problems or symptoms such as skin damage, infection, inflammation, bleeding, pigmentation may be aggravated.
- Do not use serums that contain alcohol ingredients. It may cause skin problems.
- Do not use on pets.
- Do not use near flammable substances such as thinner, benzine, polish removers or sprays.
- Do not use with anything other than the dedicated adapter and USB cord. Do not connect the dedicated adapter and USB cord to other equipment. "The cause of conflagration, electrification, or breakdown."
- Do not damage the USB cord.
- ※ Do not damage, bend more than necessary, pull, twist, place heavy objects on top of, interpose. "Causes of fire, electric shock or damage."
- Do not adopt a star-burst connection. "Causes heat generation or fire."

Prohibition of use
in a wet area

- Do not use and keep in the bathroom or high humidity places.
- Do not leave in places where water is easily splashed "on a sink etc."

Prohibition of
getting wet

- Do not immerse in water, pour or spill water (the causes of electric shock, fire or damage).

Prohibition of use
with wet hands

- Do not touch with wet hands. (the causes of electric shock or injury)

Unplug

- Make sure to hold the end of the plug when putting in or pulling the adapter. Do not grab the cord to put in or pull the adapter.
- Do not use when the adapter is damaged or the insertion to any outlet is loose.

Compulsion

- Make sure to use only household outlets.
- Insert firmly up to the root of the adapter and the plug.
- Be sure to observe the contents of the instruction manual. Failure to do so may result in some problems.

PRECAUTIONS IN USE

WARNING

Unplug

- Stop using the product, turn off the power, and pull the adapter from the outlet in the following cases.
 - Earthquakes or lightning.
 - Smoke, abnormal heat, odor or distortion of the main body, the adapter or the cord.

Prohibition

- During use, this product may slightly heat up but this is not a safety issue.
- Do not bring this product into contact with metals and any other electric appliances such as mobile phone, television.

Compulsion

- Make sure to remove all metal objects, glasses and contact lenses during use.
 - ※ Metal objects: Pierced earrings, necklaces, watches etc.
- Make sure to set to level 1 when using for the first time or until you get used to it.
- Do not press the head of product to the skin strongly and never place on the same part of the skin for longer than 3 seconds. This may cause skin troubles.
- During treatments around eyes, use this product with your eyes closed.
- Discontinue use immediately if abnormalities appear on the skin or in the body during use and consult with a doctor.

27

CAUTION

Precautions of the use of product.

Read the following sections carefully before use.

1

- Perform a special treatment to each part once a week.
 - Perform a treatment twice a week in the first and second weeks. From the third week, perform a treatment about once a week on your skin conditions.
 - Excessive treatment is not good for the skin.
- Give daily treatments to each part every day or every other day.
 - Use this product according to your skin conditions.

2

- As to facial treatments, remove makeup and wash your face to remove impurities such as makeup cosmetics and sebum on your skin and apply the dedicated beauty essence or gel, then use this product.
- Body treatment should be applied after bathing or in other clean skin conditions.

3

- Skin conditions are different depending on your physical conditions.
 - Use this product while observing your skin conditions, and discontinue use if abnormalities such as rashes, redness or itchiness appear on skin and consult with dermatologist.

PRECAUTIONS OF MAINTENANCE AND INSPECTION

⚠ WARNING

Prohibition of disassembly

- Do not disassemble, modify, or try to repair the device yourself. For repair, please inquire at dealer.

Unplug

- Unplug the adapter from the outlet when cleaning the main body. (Failure to comply may cause fire or accident.)
- Make sure to turn off the power, and disconnect the AC adapter with the outlet when not using. (Failure to comply may cause fire or damage.)

Prohibition

- Do not adhere dust and pin on AC adapter and plug. (Cause fire, electric shock or damage).
- Do not drop or hit this product.
- Do not wrap the cord around the adapter when storing the device.
- Please keep out of reach of children and do not let them use this product.
- Do not store this product in the following places.
 - Places subject to direct sunlight, high humidity and dust.
 - Places subject to high temperature or temperature below freezing point or inside a car under the blazing sun.
 - Places subject to danger of ignition, rain or wind.
 - Places subject to chemicals, gas generation or excessive vibration.
 - Places being close to heating appliances, appliances generating heat and magnetic sources.

28

Compulsion

- How to clean
Wipe this product using a soft cloth dampened with water. If wiping with a hard cloth, the product may be damaged or this may cause deterioration in performance.
- ※Do not use solvents such as alkaline agents, alcohol, thinner or benzene, polish removers, polishing powders, hot water or chemicals.
- ※Do not leave the device with the essence or gel on it after use.

INSPECTION

When this product is unused for longer than 6 months and is used again, make sure to inspect before use if there are any abnormalities (the lighting state of the LED lights and the presence of absence of cracks).

If used under faulty conditions, this may cause an accident.

NAME OF PARTS

Main Body

29

Accessories

- AC adapter
- Stand
- Instruction manual with warranty card

24KEp elements

- The surface is surrounded with 24-karat gold which is more gentle on the skin.
- ※ This metal is really soft and easily scratched. Please handle it with sufficient care. If there are any small scratches, it is not any problem with quality of the device. It is safe to use.

EXPLANATION OF FUNCTIONS

■ Special Care

In addition to TP and HP, high frequency outputs simultaneously. High-qualified home care treatment mode which pursues pleasant stimulation and the effect of warming up skin.

■ Daily Care

Enhance the efficiency of EMS by working muscles and warming skin up simultaneously.

■ DF (Duo Frequency)

2 types of high frequency are adopted.

They warm skin up from the surface to the depth of skin and lead to elastic skin.

■ TP (Tapping Pulse)

Outputting pulse intermittently realizes esthetician's hand tapping massage and pleasant stimulation by our proprietary technology. TP mode approaches various skin troubles.

■ HP (Hand Press Pulse)

Output insert pulse[※], EMS, mid-high frequency at the same time.

※Outputting insert pulse without giving any pain permeates beauty ingredients applied on skin into the skin and maximizes the effect of beauty ingredients.

■ LED

Light waves support to improve skin troubles.

■ Setting

It is possible to set volume and language, and confirm manual/ device information.

30

HOW TO USE 1. Basic operations

Screen

Bluetooth®

While connected : Brown

While unconnected : White

Start connecting with Bluetooth by tapping on Bluetooth icon.

Mode

The selected mode is displayed.

Return to previous screen

Return to previous screen by tapping on arrow icon.

31

Switch ON/OFF

■ Switch ON

Plug the AC adapter in the main body and the outlet.

■ Switch OFF

After use, unplug the AC adapter from the main body and the outlet, and then turn off the power.

Dedicated smartphone app

■ It is possible to check product information, usage history etc. with the dedicated smartphone app via Bluetooth.

※ The product can be used without being connected to the dedicated smartphone app.

◀ Please scan QR code with smartphone and download the dedicated app.

Please fill out the information according to instructions of the dedicated smartphone app.

HOW TO USE 2. Connect to dedicated app

■ Please connect Zeus III to the dedicated smartphone app.

After Turning on the power, follow the instructions below and operate Zeus III and the dedicated smartphone app in parallel.

※ If this product is not connected to the dedicated smartphone app, the usage history is not recorded.


▼ Screen of Zeus III

1 Tap on Bluetooth icon.

2 Tap on the screen.

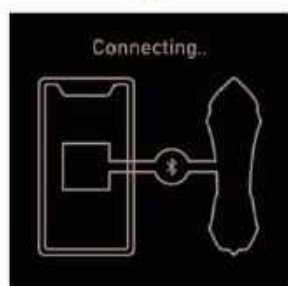▼ Screen of the dedicated smartphone app

1 Launch the dedicated smartphone app.

2 Tap on <Connect your BEAUTY DEVICE >.

3 Tap on <ACG_Zeus III >.

4 Completion screen is displayed on Zeus III and the dedicated smartphone app if the connection is completed.

※ When connecting for the first time, tap on completion screen and then register this product according to instruction of the dedicated smartphone app.

HOW TO USE 3. Setting

- 1] Swipe screen left or right and tap on Settings icon.

Volume Setting

- 1] Tap on <Volume>.

- 2] Adjust the volume by pressing Up/Down button.

33

Help

- 1] Tap on <Help>.

HOW TO USE 3. Setting

Language Setting

1 Tap on <language>.

2 Select language.

34

Confirmation of device information

1 Tap on <Device>.

It is possible to check device information such as serial number.

HOW TO USE 4. LED

1] Swipe screen left or right and tap on <LED>.

2] Select LED mode by Taping on ▲/▼.
LED mode can be used solely.

Mode 1
All colors

Mode 2
Red / Blue

Mode 3
Green / Blue

Mode 4
Flushing 3 colors

35

- Red It leads shiny and elastic skin.
- Green It leads clear and impressive skin.
- Blue It leads healthy skin.
- Lighting 3 colors For the person who has various skin troubles.
- Flushing 3 colors As the finish of treatment.

※ Do not look directly at the LED light.

※ Keep your eyes closed during the treatment around eyes.

HOW TO USE 5. Treatment

- 1 Swipe the screen left or right and select the mode from Special care, Daily Care, DF, TP, and HP by tapping on the screen.

- 2 Set power and timer
 - LED•DF(only Daily care)
 - Tap on LED to set on or off. 
 - ※ Please refer to p.34 if you would like to use LED solely.
 - Power
 - Set the power by pressing Up button and Down button.
 - ※ There are 10 power levels.
 - ※ Please use with power level 1 for the first use, and then increase power level which you feel comfortable.
 - Timer
 - Set timer by tapping on ▲/▼.
 - ※ Set timer from 1 to 10 minutes.


- 3 The selected mode starts by pressing Start/Stop button.
 - ※ A sound is made once every 1 minute.
 - ※ While mode is being used, screen can't be operated. Please press Start/Stop button when switching the mode.
- 4 Stop the mode being used by pressing Start/Stop button.
 - ※ The mode will turns off automatically when the timer runs down to 0. (The sound can be heard before the stop)

Before treatment

- 1 Remove accessories such as necklaces or rings and all metal objects that you are wearing.
- 2 Before performing a treatment, remove makeup and wash your face to remove impurities such as makeup cosmetics and sebum on your skin.
Body treatment should be applied after bathing or in other clean skin conditions.

While performing a treatment

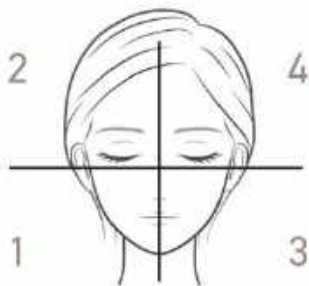- When performing a treatment for full face, do not apply the beauty essence to the face all at once but divide it into four parts to apply as illustrated left and then perform a treatment.

Daily treatment

- 《 Daily Care 》 = 10 minutes
 - Apply the gel on the treatment part and perform treatment.
- ※ Please read the treatment procedure on the page 37.

37

Special treatment once a week

- 《 Special Care 》 or 《 DF 》 or 《 TP + HP 》 = 10 minutes
- ※ Please use <Special Care>, <DF> or <TP+HP> once a week.
 - Special Care / TP / HP : Apply the beauty essence on the treatment part and perform treatment.
 - DF : Apply the gel on the treatment part and perform treatment.
- ※ Performing treatment with TP mode for 5 minutes and then performing treatment with HP mode for 5 minutes is effective.
- ※ It is possible to perform 10-minute treatment with only TP or HP.
- ※ Please read the treatment procedure on the page 38.
- ※ Treatment time is sample.


After treatment

- 1 After treatment, gently spread the beauty essence and wipe away the gel left on the skin.
- 2 After treatment, wipe away the product with a soft cloth or tissues.

TREATMENT PROCEDURE 2. Face

Daily Care

- Apply the gel on the treatment part and perform treatment.
- Place this head to the skin in order for approx. 3 seconds as shown in the illustration.
- ①~⑤ Perform each procedure two or three times.

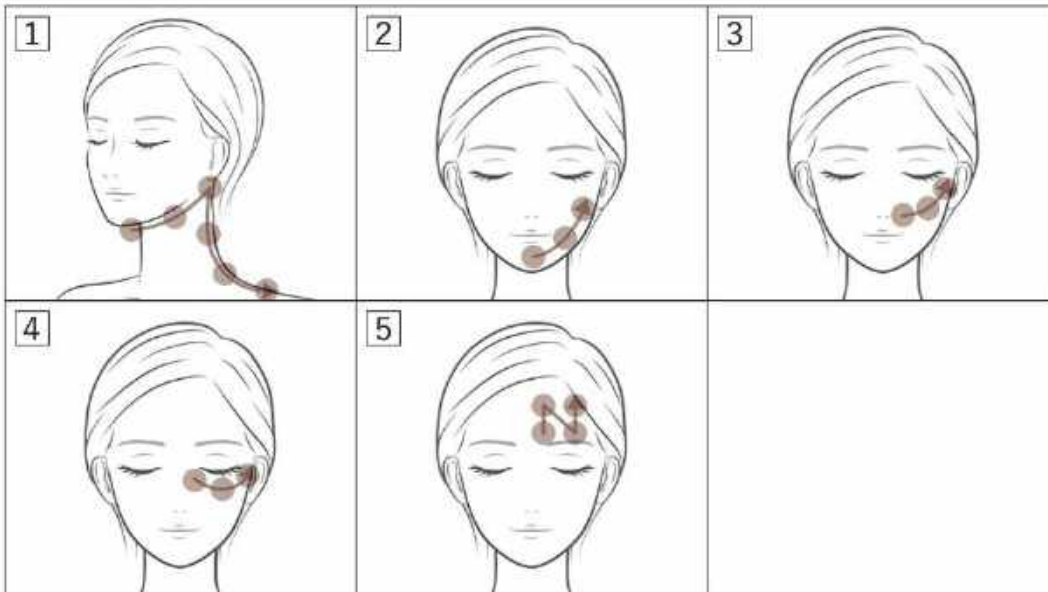- ※ The treatment procedure in the illustration is an indication.
- ※ During the treatment, apply more gel if the treatment area becomes dry. When the treatment area becomes dry, it may cause irritation on skin and the reduction of its effect.
- ※ Please adhere closely all elements on skin while treatments.
- ※ During use, wipe the surface of the product's head using a soft cloth in case the surface of the product's head becomes rough by the adhesion of dirt etc.
- ※ Do not use the product if you are physically unwell.
- ※ Excessive treatment burdens your skin. Do not use multiple times a day on the same area.


38

WARNING

Do not use on mucosa, upper eyelids, eyeballs, ears, and Adam's apple.

Special Care / DF / TP / HP


- Apply the beauty essence or gel on the treatment part and perform treatment.
- Move constantly and slowly and do not place on the same part of the skin for longer than 3 seconds.
- ①~③ Perform each procedure two or three times.

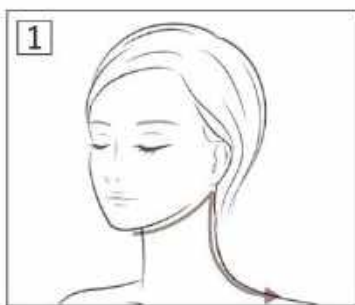Move along the face line from below the jaw to the rear side of the ear. Then, move from the parotid gland to the shoulder.

Move along the face line from the end of the jaw to the temple. Next, move along the cheek from the nostril to the temple. At the end, move along below the eye from the inner corner of the eye to the temple.

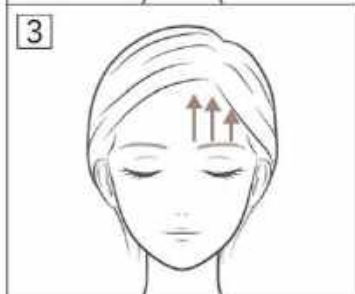Move from the bottom to the top while lifting the eyebrows upward to the hairline.

39


- ※ The treatment procedure in the illustration is an indication.
- ※ During the treatment, apply more gel or the beauty essence if the treatment area becomes dry.
When the treatment area becomes dry, it may cause irritation on skin and the reduction of its effect.
- ※ Please adhere closely all elements on skin while treatments.
- ※ During use, in the case that the surface of the product's head becomes rough by the adhesion of dirt etc., wipe the surface of the head using a soft cloth.
- ※ Do not use the product if you are physically unwell.

WARNING

Do not use on mucosa, upper eyelids, eyeballs, ears, and Adam's apple.

Upper arm


▪ Special treatment: each arm 15 minutes ▪ Daily treatment: each arm 15 minutes

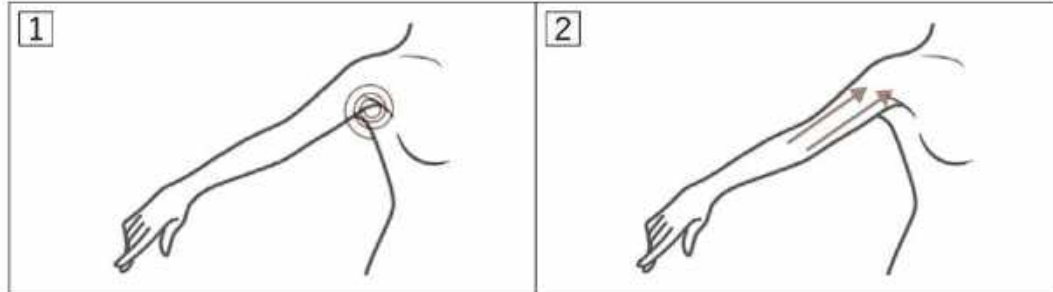- 1 Move like drawing small circles around the root of arm two or three rounds.
- 2 Move straight along the upper arm from above the elbow to the root of arm.

Waist

▪ Special treatment: 15 minutes ▪ Daily treatment: 15 minutes

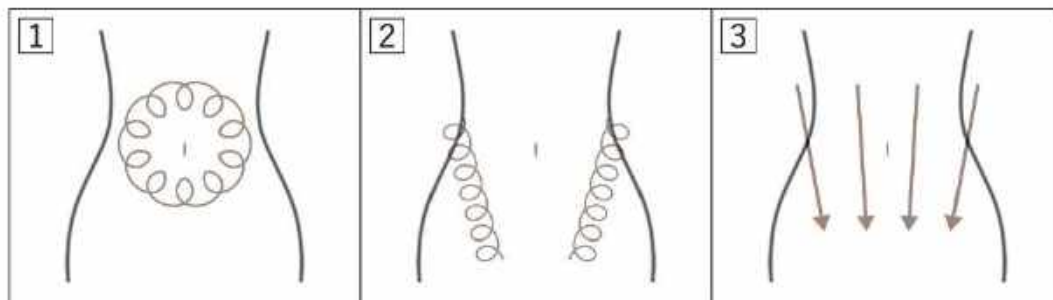
- 1 Move like drawing small circles around the navel once round.
- 2 Move like drawing small circles from the side of the waist to the inguinal node as if making the curve of the waist.
- 3 Move straight as if flowing from the top to the bottom.

40

Hip


▪ Special treatment: 15 minutes ▪ Daily treatment: 15 minutes

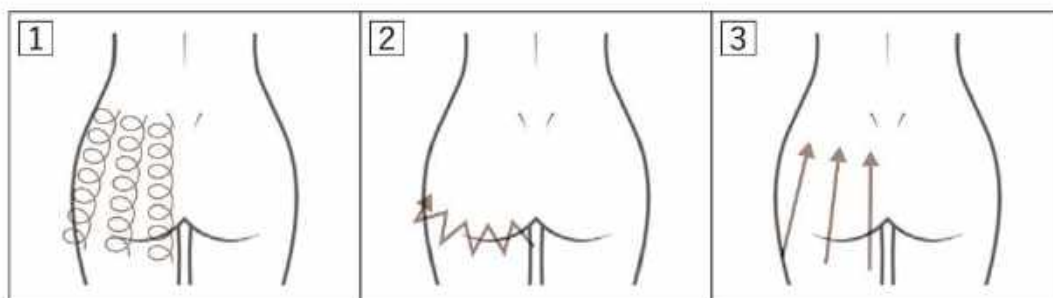- 1 Move like drawing small circles from the boundary between the hipline and the thighs to below the waist.
- 2 Move in zigzags around the boundary between the hipline and the thighs.
- 3 Move straight from the boundary between the hipline and the thighs to below the waist.

Leg (front and back)

▪ Special treatment: each calf 15 minutes ▪ Daily treatment: each calf 15 minutes

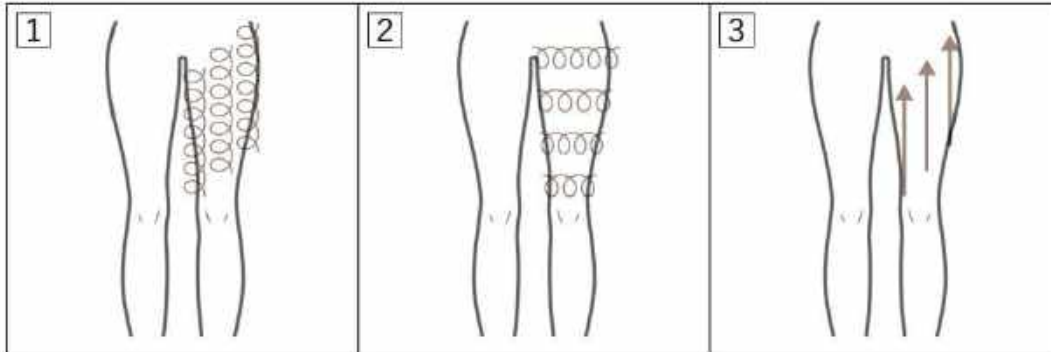
- 1 Move sideways like drawing small circles from above the knee (the back of the knee) to the root of the thighs (the hipline).
- 2 Move linearly like drawing small circles from above the knee (the back of the knee) to the root of the thighs (the hipline).
- 3 Move straight from above the knee (the back of the knee) to the root of the thighs (the hipline).

※ Perform in the same way for the front and the back of the legs.

Calf

▪ Special treatment: each calf 15 minutes ▪ Daily treatment: each calf 15 minutes

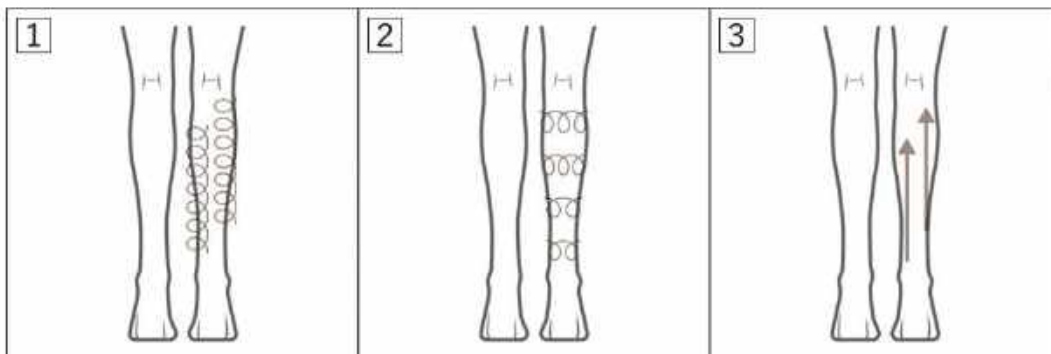- 1 Move like drawing small circles from the ankle to the back of the knee.
- 2 Move sideways like drawing small circles from the ankle to the back of the knee.
- 3 Move linearly from the ankle to the back of the knee.

PRODUCT SPECIFICATIONS

Product name	Zeus III
Power rating of AC adapter	Input 100-240V ~ 50/60Hz Output 12V = 3.0A
Power rating of Main body	Input 12V = 3.0A
Power consumption	Maximum 20W
Size	Main body: Approx. H205×W60×D66mm
Weight	Main body: Approx. 186g

42

If you have any questions, please feel free to ask.

- The way of use.
- The purchase of consumables.
- Parts exchange due to the damage, etc.
- Malfunction of the product.
- The repair of the product.

Customer Center

ARTISTIC & CO. GLOBAL Co., Ltd.

☎ 058-325-9788

Reception time: 10:00 am ~5:00 pm JST
(except Sat., Sun., and Japanese national holiday)

ASSURANCE PROVISION

■ This warranty promises free repair subject to the content written in the instruction manual.

If the product is damaged within the warranty period from the date of purchase, please contact the purchase shop or our customer center by submitting a copy of this warranty.

- ※ In the case that the serial number, date of purchase, name of dealer, customer name, address and phone number are not filled out on the warranty card, the warranty is not valid so please make sure to check that you have not missed out any required information in advance.
- ※ Attachments (AC adapter, stand, instruction manual with warranty) are not covered by the warranty.
- ※ In the case that the serial No. stickers on the main body and the warranty card are removed, the product is not covered by the warranty.

■ **There will be a charge for repairing in the following cases even within the warranty period.**

- A. In the case of the removal of the serial No. stickers of the main body and the warranty card.
- B. In the case the serial No, date of purchase, name of dealer, customer name, address and phone number are not filled out on the warranty card.
- C. Malfunction or damage due to improper handling (e.g. dropping and submerging in water) caused during transportation or moving by the user.
- D. Malfunction or damage due to fire, other natural disasters or abnormal voltage.
- E. Damage due to other devices connected to this product.
- F. Damage due to the use of an AC adapter other than the be dedicated AC adapter included with this product.
- G. In the case of repair, adjustment or modification by anyone other than us.
- H. In the case of not submitting of the warranty card.
- I. Damage due to using beauty essence or gel other than dedicated essence or gel.

■ **This warranty card is valid only in Japan.**

- ※ In the case of purchasing this product overseas, the warranty covers the product only when the store or dealer name where user purchased the product has an authorized distributor certificate issued by manufacturer.

■ **Shipping fee**

- a. The shipping fee for the return and redelivery shall be borne by our company within the warranty period. "Valid in only Japan"
- b. The shipping fee for the return and redelivery shall be borne by the customer outside the warranty period.

■ **This warranty card is not reissued. Please keep it carefully.**

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※Bluetooth® wordmark and logo are registered trademark, and ownership rights belongs to Bluetooth SIG, Inc.
※ARTISTIC&CO. GLOBAL Co., Ltd. uses these wordmark and logo under licence. Other trademark and registered trademark are every owner's.